



Topic: What is the difference between DCI 6-pin and 5-hole tubing

The **DCI** 5-hole and 6-pin use the same tubing. Both have a 6 hole terminal. The only difference is the length of the nut.

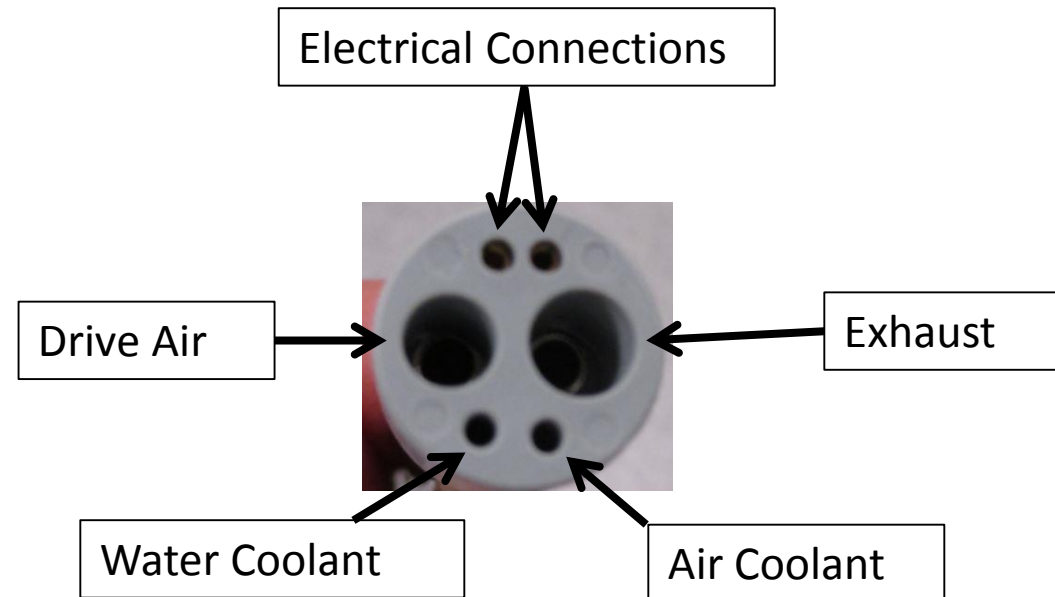

A good example that would use a 6-pin is the Kavo Lux coupler.

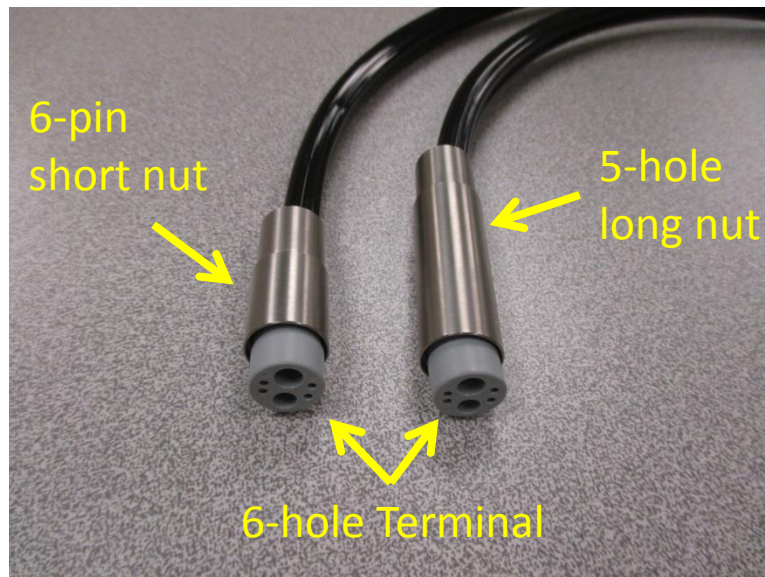

DCI Bulb
P/N 9362



6-pins



To 6-holes



The 6-pin tubing accepts a handpiece with the light bulb located in the handpiece or coupler.

Short nut



DCI 5-Hole tubing has the longer nut to accommodate the lamp module. This is used when the handpiece requires the bulb in the terminal end of the tubing.

